

SARS-CoV-2 Variants Total IgG 6-plex Panel

(Flow Cytometry Multiplex Bead Assay, Non-premixed)

【Product Name】

SARS-CoV-2 Variants Total IgG 6-plex Panel (Flow Cytometry Multiplex Bead Assay, Non-premixed)

【Catalog】

FCM-B08R

【Size】

96 Tests

【Background】

Since December 2019, the emergence of severe acute respiratory syndrome coronavirus 2 (SARS-CoV-2) and its associated disease, COVID-19, has caused a devastating pandemic worldwide. SARS-CoV-2 is a single stranded RNA virus, composed of four structural proteins, that are Spike (S) protein, Envelop (E) protein, Nucleocapsid (N) protein, and Membrane (M) protein. The viral Spike protein binds the host receptor angiotensin-converting enzyme 2 (ACE2) via the receptor binding domain (RBD), which is believed to be the mechanism that SARS-CoV-2 invades mammalian cells. SARS-CoV-2 continuously evolves, triggering pandemic waves one after another in the world. Several variants have been defined as Variant of Concern (VOC), evoking increased transmissibility and disease severity, for instance, Alpha (B.1.1.7), Beta (B.1.351), Gamma (P.1), Delta (B.1.617.2), Omicron (B.1.1.529).

Flow cytometry is originally an instrument to analyze single-cell properties in cell suspension. Nowadays, the application of flow cytometry has been extended to bead-based binding assays. Flow Cytometry Multiplex Bead Assay allows for rapid and sensitive profiling of several analytes in high-throughput format. This approach assesses dozens of analytes in a single reaction, which can generate more data with less specimen, reduced labor, and lower cost per analyte. Flow Cytometry Multiplex Bead Assay provided by ACROBiosystems are vigorously validated to ensure the test results are precise and accurate, performing robustly in supporting drug development and other intended purpose.

In principle, magnetic or non-magnetic beads with diameter 4~6 μm were colored with infrared or red fluorescent dyes, differentiated regionally by size, dye intensity, or fluorescent spectrum on flow cytometer. Each fluorescent microsphere population is conjugated to a specific target analyte, then mixed according to panel design. Processing an immunoassay in bead suspension, binding signal from bead events are acquired and statistically analyzed by flow cytometry. A common flow cytometer equipped

with blue and red lasers are sufficient to carry out multiplex immunoassay, not requiring installation of delicate instruments.

The kit, **ACROBiosystems cat#FCM-B08R**, is designed for simultaneous detection of IgG antibodies to six SARS-CoV-2 variants (Wild-type, Alpha | B.1.1.7, Beta | B.1.351, Gamma | P.1, Delta | B.1.617.2, Omicron | B.1.1.529), in a single human serum or plasma specimen.

Non-premixed beads components allow flexible combination of SARS-CoV-2 analytes.


This product is for research use only, not for use in diagnostic or therapeutic procedures.

This product provides sufficient beads suspension to carry out 96 tests, in format of 2-plex, 3-plex, 4-plex, 5-plex or 6-plex assays. **If single-plex assay is applied, only 48 tests for each analyte.**

【Principle】

SARS-CoV-2 Variants Total IgG 6-plex Panel is a multiplex bead-based indirect immunofluorescence assay, utilizing flow cytometry for data acquisition and classification. Six SARS-CoV-2 Spike RBD protein fragments (WT / Alpha / Beta / Gamma / Delta / Omicron) were conjugated to 6 fluorochrome-encoded microspheres respectively, as indicated in Table 1. The Anti-SARS-CoV-2 IgG Antibody in blood specimens was bound to the RBD encapsulated 6-plex beads, and then detected by PE-conjugated anti-human IgG detection antibody, forming bead-analyte-detection antibody sandwiches. The intensity of PE fluorescence was assessed by flow cytometry at wavelength of 575 nm approximately, in proportion to the Anti-SARS-CoV-2 IgG titer in blood specimens. The intensity of APC fluorescence was applied to classify bead populations, at wavelength of 670 nm approximately.

Table 1. Analytes and Beads

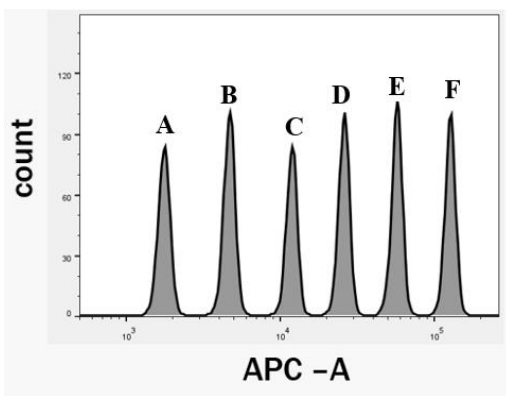

Beads	Variant of Concern
A	WT (Wild-type)
B	Alpha B.1.1.7
C	Beta B.1.351
D	Gamma P.1
E	Delta B.1.617.2
F	Omicron B.1.1.529

【Components】

Table 2. Material provided and storage condition

ID	Components	Size (96 tests)	Format	Storage (Unopened)	Storage (Opened)
FCM08-C01	SARS-CoV-2 Antibody Calibrator	20 µg	Powder	2-8°C	-20°C
FCM08-C02A	RBD WT Beads	200 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C02B	RBD Alpha Beads	200 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C02C	RBD Beta Beads	200 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C02D	RBD Gamma Beads	200 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C02E	RBD Delta Beads	200 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C02F	RBD Omicron Beads	200 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C03	Assay Buffer	40 mL×2	Liquid	2-8°C	2-8°C
FCM08-C04	5 × Wash Buffer	10 mL	Liquid	2-8°C	2-8°C
FCM08-C05	PE-Labeled Anti-human IgG Detection Antibody	200 µL	Liquid	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C06	5 × PBS Buffer	10 mL	Liquid	2-8°C	2-8°C
FCM08-C07	APC Beads	500 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C08	PE Beads	500 µL	Beads suspension	2-8°C. Light-sensitive	2-8°C. Light-sensitive
FCM08-C09	96-Well V-bottom Plate	1 plate	/	2-8°C or RT	2-8°C or RT
FCM08-C10	96-Well Sealing Film	2 pieces	/	2-8°C or RT	2-8°C or RT

【Storage and Expiration】

Unopened: Store at 2°C to 8°C upon receipt. The expiration date is labeled on the box.

Opened: The opened kit should be stored per components table. The shelf life is 30 days from the date of opening. Do not use reagents beyond expiration date.

Note: Freeze and thaw NO MORE THAN 2 times, once calibrator (ID# FCM08-C01) is reconstituted.

【Unsupplied Material and Instrument】

1. Single-channel pipettes, multi-channel pipettes, and pipette tips
2. Reagent reservoirs for multichannel pipette
3. Polypropylene microcentrifuge tubes for samples collection or dilution
4. Deionized or distilled ultrapure water
5. Bench-top centrifuge with 96-well plate swing-bucket

Optional: Vacuum manifold for 96-well plates is required, if using 96-Well Filter Plate (ACROBiosystems, cat#FCM-B01M)

6. Horizontal orbital shaker for 96-well plate
7. Vortex mixer
8. Flow cytometer equipped with two lasers: (1) Excitation at 488 nm, or 532 nm, emission around 575 nm; (2) Excitation around 633 nm, emission around 670 nm

【Important】

1. For research use only, not for use in diagnostic or therapeutic procedures.
2. Please follow the instructions strictly, for optimal and consistent data output.
3. Protect beads suspension and detection antibody from light all times to prevent photobleaching.
4. DON'T mix or substitute reagents from different kit lots. DON'T mix or substitute reagents from different manufacturers.
5. Bring the kit components to room temperature before use. Be sure that the crystal precipitates are all dissolved before use.
6. Prepare the buffer, reagent, calibrator, sample and all relevance, just prior to use.
7. Deionized or distilled water must be used for reagent preparation.
8. Ensure reagent reservoirs are clean.
9. In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test, use disposable plastic pipette tips.
10. Avoid long-term storage and repeated freeze-thaw cycles of reconstituted calibrator.

【Precaution】

All chemicals should be considered as potentially hazardous. It is recommended that this kit is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes, wash immediately with plenty of water. All blood components and biological materials should be handled and disposed properly, in accordance with local and national guideline.

【Procedure】

1. Preparation of sample and reagent

1.1 Serum collection and storage

- 1.1.1 Fresh blood samples are obtained from venous, standing at room temperature more than 30 minutes. After coagulation, centrifuge at 2000g, for 10 minutes, 4°C (Excessive centrifugation might lead to hemolysis). Carefully aspirate the serum layer, and avoid the contamination with blood cells.
- 1.1.2 Centrifuge the serum layer for 10 minutes, at 16000g, 4°C. Discard the precipitates, and keep supernatant as the freshly prepared blood serum. Use the serum immediately or keep the serum at -80°C freezer.
- 1.1.3 Dilute 5 µL freshly prepared or thawed serum in 45 µL Assay Buffer. Mix homogenously. Serum sample is ready for being tested.

1.2 Plasma collection and storage

- 1.2.1 Fresh blood samples are obtained from venous, adding anticoagulant sodium citrate, EDTA, or heparin. Centrifuge at 2000g, for 10 minutes, 4. Carefully aspirate the plasma layer, and avoid the contamination with blood cells.
- 1.2.2 Centrifuge the plasma layer 10 minutes, at 16000g, 4°C. Discard the precipitates, and keep supernatant as the freshly prepared blood plasma. Use the plasma immediately or keep the plasma at -80°C freezer.
- 1.2.3 Dilute 5 µL freshly prepared or thawed plasma in 45 µL Assay Buffer. Mix homogenously. Plasma sample is ready for being tested.

Note#1: *Frozen serum or plasma should be mixed thoroughly after thawing, and remove all visible debris by centrifuge. Thawed serum or plasma must be used up and avoid repeated freeze-thaw cycles.*

Note#2: *Hemolyzed, icteric and lipemic samples are not validated for use in this assay.*

1.3 Preparation of calibrator

1.3.1 Reconstitute 20 µg Calibrator powder in 480.7 µL Assay Buffer. Allow the bottle to sit at room temperature for 15 minutes to dissolve completely. This is calibrator stock.

Note: Mix or reconstitute protein reagent gently, avoid bubbles and foam.

1.3.2 In a microcentrifuge tube, add 270 µL Assay Buffer, 30 µL calibrator stock, and mix well. Labeled as Cal9.

1.3.3 In microcentrifuge tubes, add 150 µL Assay Buffer, labeled as Cal8, Cal7, Cal6, Cal5, Cal4, Cal3, Cal2 and Cal1, respectively.

1.3.4 Perform 2-fold serial dilutions from Cal9, as shown in table 3.

Note: Mix thoroughly before making the next dilution.

Table 3. Calibrator

Calibrator ID	Serial Dilution	Assay Buffer to add (µL)	Calibrator to add (µL)	Final Concentration (ng/mL)
Cal9	10	270	30 µL of Stock	4160
Cal8	2	150	150 µL of Cal9	2080
Cal7	2	150	150 µL of Cal8	1040
Cal6	2	150	150 µL of Cal7	520
Cal5	2	150	150 µL of Cal6	260
Cal4	2	150	150 µL of Cal5	130
Cal3	2	150	150 µL of Cal4	65
Cal2	2	150	150 µL of Cal3	32.5
Cal1	2	150	150 µL of Cal2	16.25
Cal0	-	150	-	0

1.4 Preparation of reagent and buffer

- 1.4.1** Beads suspension working solution: Define SARS-CoV-2 variants to be analyzed, and fetch corresponding RBD Beads from the kit (from ID# FCM08-C02A to ID# FCM08-C02F). Vortex RBD Beads suspension vigorously for 30 seconds or more. Immediately transfer required volume of RBD Beads to a microcentrifuge tube, and mix with required volume of Assay Buffer.

Note: Don't aspirate beads suspension less than 20 μ L to minimize pipetting errors. Perform a serial dilution if only a few beads needed for the assay.

To setup ONE test in ONE well in 96-well plate, add 2 μ L beads suspensions for EACH plex, then replenish with Assay Buffer to total volume of 220 μ L, as indicated in Table 4.

Table 4. Preparation of beads suspension working solution for one test

Multiplex	RBD Beads, ID# FCM08-C02*	Assay Buffer, ID# FCM08-C03
Single-plex	4 μ L	216 μ L
2-plex	2 μ L + 2 μ L	216 μ L
3-plex	2 μ L + 2 μ L + 2 μ L	214 μ L
4-plex	2 μ L + 2 μ L + 2 μ L + 2 μ L	212 μ L
5-plex	2 μ L + 2 μ L + 2 μ L + 2 μ L + 2 μ L	210 μ L
6-plex	2 μ L + 2 μ L + 2 μ L + 2 μ L + 2 μ L + 2 μ L	208 μ L

Note: This product provides sufficient beads suspension to carry out 96 tests, in format of 2-plex, 3-plex, 4-plex, 5-plex or 6-plex assays.

If single-plex assay is applied, only 48 tests for each analyte.

- 1.4.2** Detection antibody working solution: Bring Detection Antibody to room temperature. Aspirate 200 μ L Detection Antibody, mixed with 20 mL Assay Buffer.
- 1.4.3** 1 \times wash buffer: Bring 5 \times Wash Buffer to room temperature. Aspirate 10 mL 5 \times Wash Buffer, mixed with 40 mL ultrapure water.
- 1.4.4** 1 \times PBS: Bring 5 \times PBS Buffer to room temperature. Aspirate 10 mL 5 \times PBS Buffer, mixed with 40 ml ultrapure water.

2. Assay procedure

Note: Performing the assays with **96-Well V-bottom Plate (ACROBiosystems, cat#FCM-B08R, ID# FCM08-C09)**.

- 2.1** Add serial dilutions of calibrator, or properly diluted blood samples to 96-Well V-bottom Plate, 30 μ L per well.

Note: Run calibrators in duplicates. Follow the attached Plate Layout to achieve good accuracy.

- 2.2** Add beads suspension working solution to 96-well plate, 220 μ L per well.
2.3 Seal the plate. Incubate at room temperature for 60 minutes, with continuous shaking to ensure the beads always suspended homogenously in the solution. Avoid light.

Note: Wrap the 96-well plate by aluminum foil, or keep orbital shaker in the dark.

- 2.4** Centrifuge at 3000g, for 10 minutes, at room temperature. Discard supernatant.
2.5 Add 200 μ L 1 \times wash buffer to each well. Mix by pipetting up and down for 6-8 times.
2.6 Centrifuge at 3000g, for 10 minutes, at room temperature. Discard wash buffer.
2.7 Add detection antibody working solution, 200 μ L per well.
2.8 Seal the plate. Incubate at room temperature for 40 minutes, with continuous shaking to ensure the beads always suspended homogenously in the solution. Avoid light.
2.9 Centrifuge at 3000g, for 10 minutes, at room temperature. Discard supernatant.
2.10 Add 200 μ L 1 \times PBS to each well. Mix by pipetting up and down. Ensure the beads well separated and not aggregated.
2.11 Subject to flow cytometry analysis. If not being analyzed immediately, store at 2~8 $^{\circ}$ C and protect from light. Read by flow cytometry within 2 hours.

Note: Resuspend beads immediately prior to reading by pipetting up and down.

2*. Assay procedure

Note: Performing the assays with **96-Well Filter Plate (ACROBiosystems, cat#FCM-B01M)**.

- 2.1*** Add serial dilutions of calibrator, or properly diluted blood samples to 96-Well Filter Plate, 30 μ L per well.

Note: Run calibrators in duplicates. Follow the attached Plate Layout to achieve good accuracy.

- 2.2*** Add beads suspension working solution to 96-well plate, 220 μ L per well.
2.3* Seal the plate. Incubate at room temperature for 60 minutes, with continuous shaking to ensure the beads always suspended homogenously in the solution. Avoid light.

Note: Wrap the 96-well plate by aluminum foil, or keep orbital shaker in the dark.

- 2.4*** Place the 96-well plate on the manifold. Remove the sealing film and remove supernatant by vacuum.
2.5* Add 200 μ L 1 \times wash buffer to each well. Mix by pipetting up and down for 6-8 times.

- 2.6* Place the 96-well plate on the manifold and remove wash buffer by vacuum.
- 2.7* Add detection antibody working solution, 200 μ L per well.
- 2.8* Seal the plate. Incubate at room temperature for 40 minutes, with continuous shaking to ensure the beads always suspended homogenously in the solution. Avoid light.
- 2.9* Place the 96-well plate on the manifold. Remove the sealing film and remove supernatant by vacuum.
- 2.10* Add 200 μ L 1 \times PBS to each well. Mix by pipetting up and down. Ensure the beads well separated and not aggregated.
- 2.11* Subject to flow cytometry analysis. If not being analyzed immediately, store at 2~8 $^{\circ}$ C and protect from light. Read by flow cytometry within 2 hours.

Note: Resuspend beads immediately prior to reading by pipetting up and down.

3. Flow cytometer setup

- 3.1 Flow cytometer equipped with two lasers are compatible with the assay: (1) excitation laser at 488 nm or 532 nm, and emission around 575 nm; (2) excitation laser around 633 nm, and emission around 670 nm. Instruments tested by this assay were represented in table 5.

Table 5. Partial list of compatible flow cytometers

Manufacturer	Flow Cytometer	Classification Channel	Reporter Channel
BD Biosciences	BD FACSLyric™	APC	PE
Beckman Coulter	Cytoflex S	APC	PE
Thermo Fisher Scientific	Attune NxT	RL1	YL1
Luminex Corporation	Guava easyCyte3L	RED-R	YEL-B

- 3.2 Start up and calibrate flow cytometer following the instrument user guide.
- 3.3 Enable FSC and SSC channel.
- 3.4 Enable the fluorescent channel with an excitation peak at 488 nm or 532 nm, and an emission peak around 575 nm (for example, PE channel). Enable the fluorescent channel with an excitation peak at 633 nm, and an emission peak around 670 nm (for example, APC channel).
- 3.5 Load APC Beads, and adjust PMT gain or voltage of FSC and SSC, positioning the singlet beads in the central region of FSC-SSC scatter plot. Adjust FSC and SSC threshold to diminish debris.
- 3.6 In FSC-SSC scatter plot, draw a region around the single beads to exclude doublets or more, and gate all subsequent analysis on this region, labeled as P1.
- 3.7 In APC-count histogram, P1 sub-population represents two peaks. Adjust gain or voltage of

- APC channel, positioning APC positive peak close to upper limit of APC scale (around 10^5), and APC negative peak between 0 ~ 1000 correspondingly.
- 3.8** Remove APC Beads and load PE Beads. In FSC-SSC scatter plot, draw P1 gate in the same way as above. In PE-count histogram, P1 sub-population represents two peaks. Adjust gain and voltage of PE channel, positioning PE positive peak close to upper limit of PE scale (around 10^5), and PE negative peak between 0 ~ 1000 correspondingly. Thus, the gain or voltage is configured.
 - 3.9** Load APC Beads and PE Beads once again to adjust fluorescence spillover of APC and PE channel, following the instrument user guide or former experiences to establish compensation matrix. (In most instances, APC doesn't spillover into PE channel, and vice versa).
 - 3.10** Select medium flow rate.
 - 3.11** Set up stop criteria as 4800 events or beads collected in P1 gate.

Note: It is highly recommended to apply the calibrator sample of the highest concentration (data point of 4160 ng/ml) to verify the instrument and method configuration. Make further adjustments if necessary.

4. Data Acquisition and Analysis


4.1 Data acquisition

- 4.1.1** Make sure the flow cytometer is well tuned, following the instrument user guide and method configuration illustrated above.
- 4.1.2** Create an experiment in 96-well plate format.

Note: If 96-well plate loader is not available, transfer the samples to FACS tubes, replenish 100 μ L $1 \times$ PBS, and read one by one.

- 4.1.3** Resuspend beads by pipetting up and down.
- 4.1.4** Load the plate and start acquisition.
- 4.1.5** In APC-count histogram, create P2 ~P7 gates in P1 sub-population. Six-plex RBD beads was differentiated, as indicated in table 6.
- 4.1.6** Record mean fluorescence intensity (MFI) of PE channel.

Table 6. Segregation of analyte-specific populations

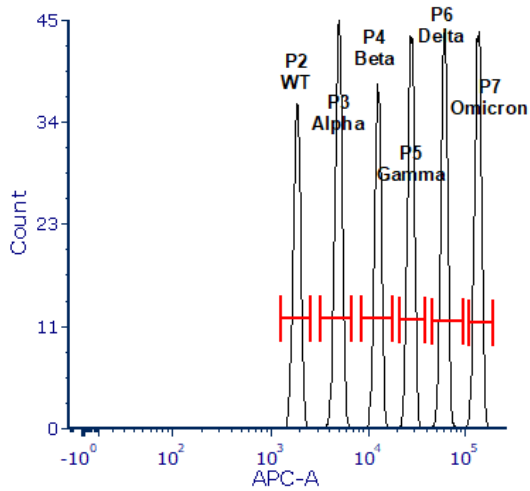


Beads	Variant of Concern	Gate
A	WT (Wild-type)	P2
B	Alpha B.1.1.7	P3
C	Beta B.1.351	P4
D	Gamma P.1	P5
E	Delta B.1.617.2	P6
F	Omicron B.1.1.529	P7

4.2 Data analysis

- 4.2.1** Construct a calibration curve for each analyte by plotting Log₁₀ concentration value of serial diluted calibrators against mean fluorescence intensity (MFI) of PE channel. A four parameter logistic (4-PL) curve model is applied, with R² value above 0.99.
- 4.2.2** Calculate the concentration of unknown from the calibrator curve of each analyte.
- 4.2.3** Determine the concentration of blood specimens after multiplying by the dilution factor (10×).

【Quick Guide】

Print the Quick Guide and follow the protocol.

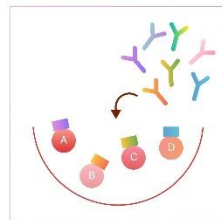

1


- ★ Add 30 μ L calibrator or blood samples
- ★ Add 220 μ L beads suspension working solution to 96-well plate

Shake, avoid light
18-25°C 60 minutes



2


- ★ Wash with 200 μ L 1 \times wash buffer, pipetting up and down for 6-8 times.

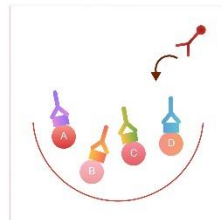

3

- ★ Add 200 μ L detection antibody working solution

Shake, avoid light
18-25°C 40 minutes

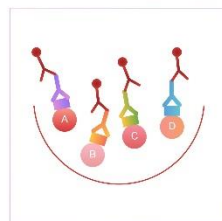

4

- ★ Resuspend with 200 μ L 1 \times PBS

avoid light <2hours



5

- ★ Flow cytometry analysis

【Troubleshooting】

Concerns	Possibilities	Suggestions
After centrifuge, the white bead precipitates are not visible or become less and less after multiple-step operation.	The pellets are very loosely attached to the well, and lost during aspiration.	Aspirate the supernatant slowly and carefully. Don't let the plate sit very long. Centrifuge again if necessary.
Variation of bead count in duplicated wells.	Aspiration takes so long time that the beads settled to the bottom of the tube.	Quickly aspirate and dispense the bead suspension. Vortex beads vigorously before first use, and vortex briefly in between.
Plenty of debris were observed in FSC-SSC scatter plot during data acquisition.	Improper setting of FSC and SSC threshold.	Increase threshold value of FSC and SSC.
Plenty of bead doublets are observed by plotting FSC height versus FCS area.	Beads aggregate due to long time sitting or insufficient resuspending.	Resuspend the beads by pipetting up and down vigorously, then re-load onto flow cytometer.
Less than 6 bead populations in APC-count histogram.	The PMT gain or voltage value of APC fluorescent channel is too high.	Adjust PMT gain or voltage of APC fluorescent channel, ensuring 6 intact peaks observed.
Less than 6 bead populations in APC-PE scatter plot, though all the 6 peaks obtained in APC-count histogram.	The PMT gain or voltage value of PE fluorescent channel is too high.	Adjust PMT gain or voltage of PE fluorescent channel, ensuring 6 bead populations in APC-PE scatter plot.
PE fluorescent intensity of low concentration calibrator is higher than that of high concentration calibrator.	Insufficient needle wash and clean between samples.	At least one washing cycle between samples in flow cytometry setting.
		Follow the Plate Layout suggested, and read the plate by columns to reduce cross-contamination on flow cytometer.

【Plate Layout】

Arrange the samples in vertical, and read the plate by column.

	1	2	3	4	5	6	7	8	9	10	11	12
A	Blank	Cal4	Cal8								
B	Blank	Cal4	Cal8								
C	Cal1	Cal5	Cal9								
D	Cal1	Cal5	Cal9								
E	Cal2	Cal6	Sample 1								
F	Cal2	Cal6	Sample 1								
G	Cal3	Cal7	Sample 2								
H	Cal3	Cal7	Sample 2								