SARS-CoV-2 Spike Trimer Protein (T19R, V70F, FR157-158Del, A222V, W258L, K417N, L452R, T478K, D614G, P681R, D950N), His Tag (MALS verified)

Catalaa # ODNI OEOLII

Synonym

Spike,S protein,Spike glycoprotein,S glycoprotein

Source

SARS-CoV-2 Spike Trimer, His Tag (SPN-C52Ht) is the ectodomain of SARS-CoV-2 spike protein which contains AA Val 16 - Pro 1213 (Accession # <u>QHD43416.1</u>). Mutations T19R, V70F, FR157-158Del, A222V, W258L, K417N, L452R, T478K, D614G, P681R, D950N were identified in the SARS-CoV-2 Delta sub-lineages (Pango lineage: AY.1/AY.2/AY.3; also known as "Delta Plus"). The recombinant protein is expressed from human 293 cells (HEK293) with T4 fibritin trimerization motif and a polyhistidine tag at the C-terminus. Proline substitutions (F817P, A892P, A899P, A942P, K986P, V987P) and alanine substitutions (R683A and R685A) are introduced to stabilize the trimeric prefusion state of SARS-CoV-2 S protein and abolish the furin cleavage site, respectively.

Predicted N-terminus: Val 16

Molecular Characterization


This protein carries a polyhistidine tag at the C-terminus

The protein has a calculated MW of 137.8 kDa. The protein migrates as 160-200 kDa under reducing (R) condition (SDS-PAGE) due to glycosylation.

Endotoxin

Less than 1.0 EU per μ g by the LAL method.

SDS-PAGE

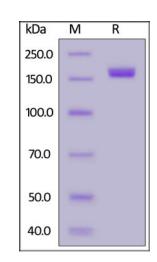

SARS-CoV-2 Spike Trimer, His Tag on SDS-PAGE under reducing (R) condition. The gel was stained overnight with Coomassie Blue. The purity of the protein is greater than 95%.

Purity

>95% as determined by SDS-PAGE.

Formulation

Lyophilized from 0.22 µm filtered solution in PBS with trehalose as protectant.

Contact us for customized product form or formulation.

Reconstitution

Please see Certificate of Analysis for specific instructions.

For best performance, we strongly recommend you to follow the reconstitution protocol provided in the CoA.


Storage

For long term storage, the product should be stored at lyophilized state at -20°C or lower.

Please avoid repeated freeze-thaw cycles.

This product is stable after storage at:

- -20°C to -70°C for 12 months in lyophilized state;
- -70°C for 3 months under sterile conditions after reconstitution.

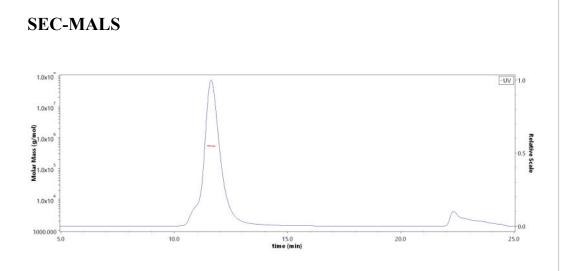

The purity of SARS-CoV-2 Spike Trimer, His Tag (Cat. No. SPN-C52Ht) is more than 85% and the molecular weight of this protein is around 510-564 kDa verified by SEC-MALS.

<u>Report</u>


Bioactivity-ELISA

>>> www.acrobiosystems.com

4/10/2023


SARS-CoV-2 Spike Trimer Protein (T19R, V70F, FR157-158Del, A222V, W258L, K417N, L452R, T478K, D614G, P681R, D950N), His Tag (MALS verified)



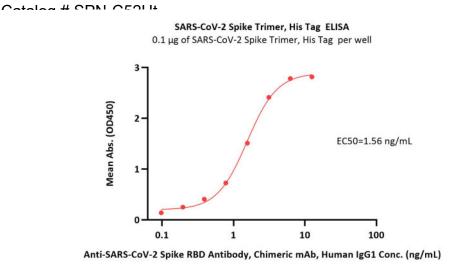


SARS-CoV-2 Spike Trimer, His Tag Conc. (ng/mL)

Immobilized SARS-CoV-2 Spike Trimer, His Tag (Cat. No. SPN-C52Ht) at 1 μ g/mL (100 μ L/well) can bind Anti-SARS-CoV-2 Spike RBD Antibody, Chimeric mAb, Human IgG1 (Cat. No. S1N-M122) with a linear range of 0.1-3 ng/mL (QC tested).

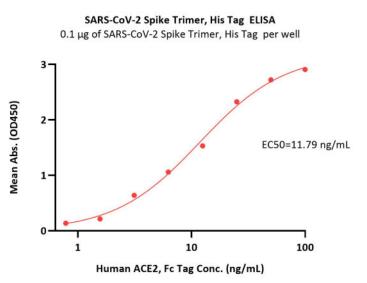

Immobilized SARS-CoV-2 Spike Trimer, His Tag (Cat. No. SPN-C52Ht) at 1 μ g/mL (100 μ L/well) can bind Human ACE2, Fc Tag (Cat. No. AC2-H5257) with a linear range of 0.8-13 ng/mL (Routinely tested).

Immobilized Anti-SARS-CoV-2 Spike RBD Antibody, Chimeric mAb, Human IgG1 (Cat. No. S1N-M122) at 1 μ g/mL (100 μ L/well) can bind SARS-CoV-2 Spike Trimer, His Tag (Cat. No. SPN-C52Ht) with a linear range of 0.1-3 ng/mL (Routinely tested).

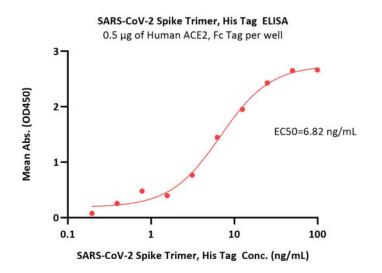

Immobilized Human ACE2, Fc Tag (Cat. No. AC2-H5257) at 5 μ g/mL (100 μ L/well) can bind SARS-CoV-2 Spike Trimer, His Tag (Cat. No. SPN-C52Ht) with a linear range of 0.2-13 ng/mL (Routinely tested).

Background

It's been reported that SARS-CoV-2 can infect the human respiratory epithelial cells through interaction with the human ACE2 receptor. The spike protein is a large type I transmembrane protein containing two subunits, S1 and S2. S1 mainly contains a receptor binding domain (RBD), which is responsible for recognizing the cell surface receptor. S2 contains basic elements needed for the membrane fusion. The S protein plays key parts in the induction of neutralizing-antibody and T-cell responses, as well as protective immunity.

Clinical and Translational Updates

Please contact us via <u>TechSupport@acrobiosystems.com</u> if you have any question on this product.



4/10/2023