

#### Synonym

Spike, Sprotein, Spike glycoprotein, Sglycoprotein

#### Source

Biotinylated SARS-CoV-2 Spike Trimer, His,Avitag<sup>TM</sup> (BA.2/Omicron) (SPN-C82Er) is expressed from human 293 cells (HEK293). It contains AA Val 16 - Pro 1213 (Accession # QHD43416.1 (T19I, LPP24-26del, A27S, G142D, V213G, G339D, S371F, S373P, S375F, T376A, D405N, R408S, K417N, N440K, S477N, T478K, E484A, Q493R, Q498R, N501Y, Y505H, D614G, H655Y, N679K, P681H, N764K, D796Y, Q954H, N969K, R683A, R685A, F817P, A892P, A899P, A942P, K986P, V987P)). The spike mutations are identified on the SARS-CoV-2 Omicron variant (Pango lineage: BA.2; GISAID clade: GRA; Nextstrain clade: 21L). The recombinant protein is expressed from human 293 cells (HEK293) with T4 fibritin trimerization motif and a polyhistidine tag at the C-terminus. Proline substitutions (F817P, A892P, A899P, A942P, K986P, V987P) and alanine substitutions (R683A and R685A) are introduced to stabilize the trimeric prefusion state of SARS-CoV-2 S protein and abolish the furin cleavage site, respectively.

Predicted N-terminus: Val 16

#### **Molecular Characterization**

This protein carries a polyhistidine tag at the C-terminus, followed by an Avi tag (Avitag<sup>TM</sup>)

The protein has a calculated MW of 139.7 kDa. The protein migrates as 160-190 kDa under reducing (R) condition (SDS-PAGE) due to glycosylation.

# Labeling

Biotinylation of this product is performed using Avitag<sup>TM</sup> technology. Briefly, the single lysine residue in the Avitag is enzymatically labeled with biotin.

### **Protein Ratio**

Passed as determined by the HABA assay / binding ELISA.

#### Endotoxin

Less than 1.0 EU per μg by the LAL method.

### **Purity**

>95% as determined by SDS-PAGE.

>90% as determined by SEC-MALS.

#### **Formulation**

Lyophilized from 0.22 µm filtered solution in PBS with trehalose as protectant.

Contact us for customized product form or formulation.

#### Reconstitution

Please see Certificate of Analysis for specific instructions.

For best performance, we strongly recommend you to follow the reconstitution protocol provided in the CoA.

### Storage


For long term storage, the product should be stored at lyophilized state at -20 $^{\circ}$ C or lower.

Please avoid repeated freeze-thaw cycles.


This product is stable after storage at:

- -20°C to -70°C for 12 months in lyophilized state;
- -70°C for 3 months under sterile conditions after reconstitution.

### **SDS-PAGE**

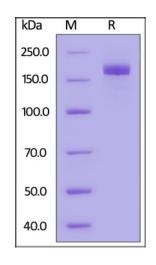

# SEC-MALS

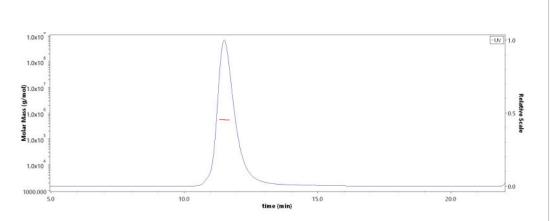

# Biotinylated SARS-CoV-2 Spike Trimer Protein, His,Avitag™ (BA.2/Omicron) (MALS verified)

Catalog # SPN-C82Er




Biotinylated SARS-CoV-2 Spike Trimer, His, Avitag (BA.2/Omicron) on SDS-PAGE under reducing (R) condition. The gel was stained overnight with Coomassie Blue. The purity of the protein is greater than 95%.

The purity of Biotinylated SARS-CoV-2 Spike Trimer, His,Avitag (BA.2/Omicron) (Cat. No. SPN-C82Er) is more than 90% and the molecular weight of this protein is around 525-575 kDa verified by SEC-MALS.

Report


### **Bioactivity-ELISA**

**Biotinylated SARS-CoV-2 Spike Trimer, His, Avitag (BA.2/Omicron) ELISA**0.1 μg of Biotinylated SARS-CoV-2 Spike Trimer, His, Avitag (BA.2/Omicron) per well

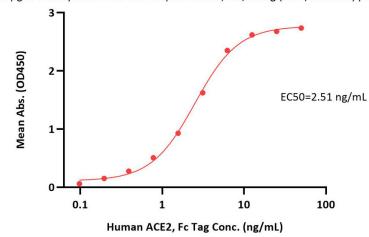

Immobilized Biotinylated SARS-CoV-2 Spike Trimer, His,Avitag (BA.2/Omicron) (Cat. No. SPN-C82Er) at 1 μg/mL (100 μL/well) on streptavidin (Cat. No. STN-N5116) precoated (0.5 μg/well) plate can bind Human ACE2, Fc Tag (Cat. No. AC2-H5257) with a linear range of 0.1-6 ng/mL (QC tested).


Biotinylated SARS-CoV-2 Spike Trimer, His,Avitag (BA.2/Omicron) ELISA 0.1  $\mu$ g of Biotinylated SARS-CoV-2 Spike Trimer, His,Avitag (BA.2/Omicron) per well

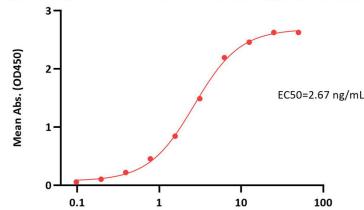

Anti-SARS-CoV-2 Spike RBD Antibody, Chimeric mAb, Human IgG1 Conc. (ng/mL)

Immobilized Biotinylated SARS-CoV-2 Spike Trimer, His,Avitag (BA.2/Omicron) (Cat. No. SPN-C82Er) at 1 μg/mL (100 μL/well) on streptavidin (Cat. No. STN-N5116) precoated (0.5 μg/well) plate can bind Anti-SARS-CoV-2 Spike RBD Antibody, Chimeric mAb, Human IgG1 (Cat. No. S1N-M122) with a linear range of 0.1-6 ng/mL (Routinely tested).

**Biotinylated SARS-CoV-2 Spike Trimer, His, Avitag (BA.2/Omicron) ELISA** 0.1 μg of Biotinylated SARS-CoV-2 Spike Trimer, His, Avitag (BA.2/Omicron) per well

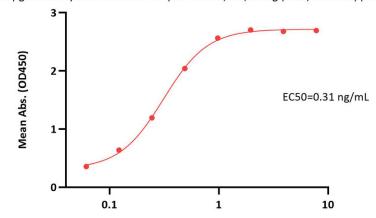

Anti-SARS-CoV-2 Spike RBD Antibody, Mouse IgG1 (Omicron Specific) Conc. (ng/mL)

Immobilized Biotinylated SARS-CoV-2 Spike Trimer, His,Avitag (BA.2/Omicron) (Cat. No. SPN-C82Er) at 1 μg/mL (100 μL/well) on streptavidin (Cat. No. STN-N5116) precoated (0.5 μg/well) plate can bind Anti-SARS-CoV-2 Spike RBD Antibody, Mouse IgG1 (Omicron Specific) (Cat. No. SPD-M305) with a linear range of 0.1-1 ng/mL (Routinely tested).

### Background

It's been reported that SARS-CoV-2 can infect the human respiratory epithelial cells through interaction with the human ACE2 receptor. The spike protein is a large type I transmembrane protein containing two subunits, S1 and S2. S1 mainly contains a receptor binding domain (RBD), which is responsible for recognizing the cell surface receptor. S2 contains basic elements needed for the membrane fusion. The S protein plays key parts in the induction of neutralizing-antibody and T-cell responses, as well as protective immunity.

## **Clinical and Translational Updates**

# Biotinylated SARS-CoV-2 Spike Trimer Protein, His,Avitag™ (BA.2/Omicron) (MALS verified)





Please contact us via <u>TechSupport@acrobiosystems.com</u> if you have any question on this product.